

# G.M. 1: Basics of GM Counter-Characteristics and Counting Statistics

Aritra Mukhopadhyay  
National Institute of Science Education and Research  
Bhubaneswar, Odisha 751005, India  
3rd year, Integrated M.Sc. Physics  
Roll No.: 2011030  
(Dated: April 6, 2023)

The experiment aimed to explore the basics of GM Counter-Characteristics and Counting Statistics in nuclear physics. The Geiger-Muller (GM) counter was used to study radioactivity, determine its characteristics, confirm the inverse square law related to gamma rays, calculate source efficiency, and analyze counting statistics. The results showed that the GM characteristic curve closely resembled an ideal one and confirmed the inverse square law. It was also found that the GM counter was more efficient in detecting radiation from  $Tl^{204}$  ( $\beta$ -rays). However, the counting statistics distributions were incomplete due to insufficient sample numbers or defects.

## I. THEORY

### A. Description of the Geiger-Muller counter

The Geiger-Muller (GM) counter works by detecting the ion-electron pairs created by the interaction of charged particles in a gas mixture. The GM tube is a metal cylinder with a thin wire (anode) at its axis and a metal cylinder (cathode) maintained at a high voltage to create an electric field. The radiation enters through a window on the tube, creating ion-electron pairs that are swept by the electric field to produce a phenomenon called an avalanche, which generates output pulses that are counted by related circuits. A schematic diagram of the GM counter is shown in Figure 1.

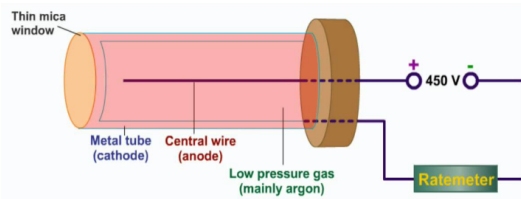


FIG. 1: Typical GM counter characteristics

As shown in Figure 2 the operating voltage of the GM counter is set in the plateau region, where the counting rate is relatively constant. The plateau length and slope determine the stability of the counting rate, while the dead time, resolving time, and recovery time limit the counting rate. The GM counter is insensitive to ionizing events during the dead time, and the resolving time and recovery time set the minimum time interval between two distinct and normal-size pulses, respectively. Higher voltages and the gas composition inside the GM tube can reduce the effects of these factors. The background counting rate is due to cosmic rays or other active sources in the same room.

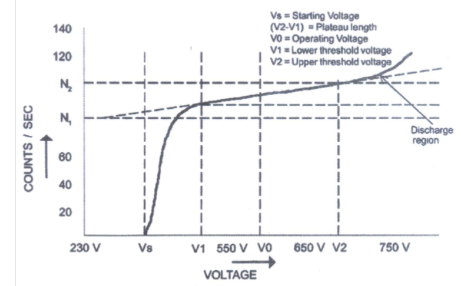


FIG. 2: Typical GM counter characteristics

### B. Inverse square law: Gamma rays

It states that the gamma radiations reduce inversely proportional to the distance,  $D$ , between the detector and the radiation source. Thus, the counting rate,  $R$  (counts/second), should be related as follows:

$$R \propto \frac{1}{D^2}$$

$$\therefore \log(R) = -2 \times \log(D) + c$$

where  $c$  is some constant

### C. Efficiency of GM counter

Knowing the activity of that specific source allows one to calculate the effectiveness of a GM counter. The quantity of the radiation source's disintegrations per second is known as activity. The ratio of the observed counts per second to the number of disintegrations that are detected by the detector each second is referred to as the efficiency. The efficiency,  $E$ , is now:

$$E = \frac{CPS}{DPS} \times 100 \quad (1)$$

where CPS is the counts per second, and DPS is the disintegrations per second, falling on the window of the detector. DPS can be calculated as:

$$DPS = A \frac{d^2}{16D^2} \quad (2)$$

where A is the activity of the source, d is the diameter of the detector, and D is the distance between the source and the detector.

#### D. Counting Statistics

Say  $N_i$  denote the  $i^{th}$  reading of a measurement in a set of  $n$  measurements, then the equations for calculating mean ( $\bar{N}$ ), variance  $\sigma^2$ , and standard deviation  $\sigma$  (for large samples) are:

$$\bar{N} = \frac{1}{n} \sum_{i=1}^n N_i$$

$$\sigma^2 = \frac{1}{n} \sum_{i=1}^n (\bar{N} - N_i)^2$$

## II. OBSERVATIONS & CALCULATIONS

#### A. GM characteristics

A gamma source ( $^{137}\text{Cs}$ ) was used and the voltage was varied to get the data in Table 1. From the data, the curve was plotted as shown in Figure 3.

Potential (V)	Count	Background	Corrected Counts
343	0	0	0
344	2342	0	2342
360	4302	34	4268
390	4963	53	4910
420	4939	53	4886
450	5114	33	5081
480	5081	38	5043
510	5327	40	5287
540	5264	44	5220
570	5264	36	5228
600	5400	45	5355
630	5394	57	5337
660	9570	87	9483
690	10320	76	10244

TABLE I: GM characteristic data

From the graph we can see:

1. Lower threshold volatge = 390V
2. Upper threshold voltage = 630V

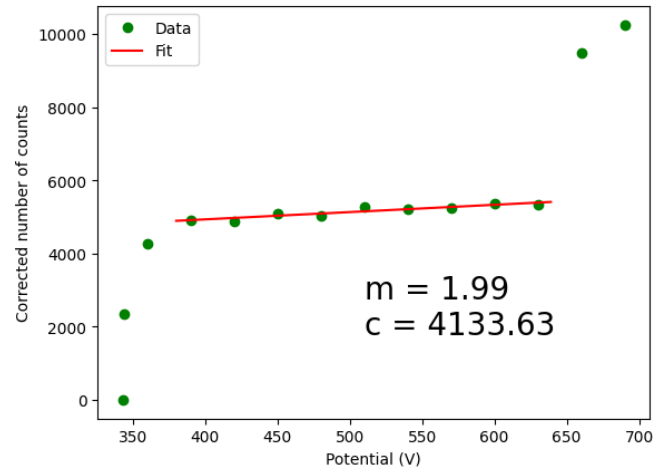


FIG. 3: Graph of GM characteristics

3. Plateu length =  $630 - 390 = 240V$
4. Counts at lower threshold voltage = 4910
5. Slope = 1.99
6. Plateau Slope percent ( $s$ ) =  $\frac{\text{slope}}{N} * 10000 = \frac{1.99}{4910} * 10000 = 4.073\%$
7. Operating Volatage =  $\frac{600+360}{2} = 510V$

#### B. Inverse Square Law

Background	Average
576	558
582	
600	
486	
564	
540	

TABLE II: Background Count (60s)

Distance (d) (cm)	Count	Corrected Count in 60s	Net Count Rate (R) (per second)	Product $C = Rd^2$	$\log(d)$	$\log(R)$
2.0	10462	9904	165.067	660.27	0.693	5.106
2.5	7797	7239	120.650	754.06	0.916	4.792
3.0	5819	5261	87.683	789.15	1.098	4.473
3.5	4515	3957	65.950	807.89	1.252	4.188
4.0	3624	3066	51.100	817.60	1.386	3.933
4.5	2961	2403	40.050	811.01	1.504	3.690
5.0	2424	1866	31.100	777.50	1.609	3.437
5.5	2062	1504	25.067	758.27	1.704	3.221
6.0	1824	1266	21.100	759.60	1.791	3.049
7.0	1470	912	15.200	744.80	1.945	2.721

TABLE III: Inverse Square Law Data & Calculations

[illegible]

